


The Thermo Scientific VisionMate Wireless 2D Barcode Reader decodes both 2D barcodes, as well as linear barcodes. Use it as a wireless reader via Bluetooth® technology for lab, cold room or in-field sample reading onto your laptop, PC or PDA.

Thermo Scientific VisionMate® Wireless 2D Barcode Reader

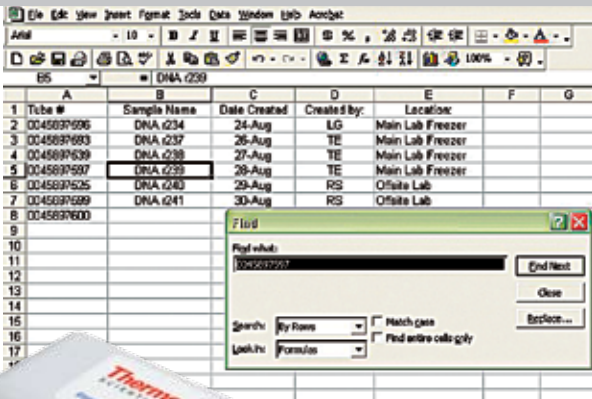

Scanning Simplicity and Flexibility.

Hold a Thermo Scientific Matrix® 2D Barcoded Storage Tube over the Thermo Scientific Matrix VisionMate Wireless' scanning window and the code is decoded and imported instantly into any application where you can place your cursor, including spreadsheets, databases and LIMS. Just as easily scan and record linear (1D) barcodes into your applications. Also enable barcode reading in your cold room, in the field or across the laboratory via Bluetooth technology.

Effortless Operation. The instrument's flat bottom profile allows for scanning without the need to hold the reader, yet at only 13.5 ounces (383 grams) the VisionMate Wireless is easy to handle.

Knowledgeable Support. Our team has a strong reputation for supporting both equipment and consumables, and can demonstrate simple ways to track samples in your laboratory using our VisionMate Wireless reader.

Thermo Scientific VisionMate Wireless 2D Barcode Reader

Dimensions (W x D x H)	1.68"/18.8 cm x 7.4"/4.3 cm x 2.28"/5.8 cm
Weight	13.5 oz/383 g
Battery Type	1300 mA Lithium Ion Battery
Battery Run Time	> 8 hours of continuous run time
Charging Time	Approximately 6 hours
Charging Connection	Via supplied USB port connection
Battery Indicator Light	Green: 100-50% charged Amber: 50-20% charged Red: 20-0% charged
Supported Barcode Types	Any Thermo Scientific Matrix 2D Barcoded Storage Tube and virtually any other 2D barcode, as well as any linear (1D) code
Computer Requirements	PC/Laptop: Microsoft® Windows® 2000 or Windows XP with available USB 1.1 port PDA: Pocket PC w/ Active-Sync (built-in Bluetooth required)
Connection Methods	Wireless: via supplied USB Bluetooth Adapter or via built-in Bluetooth Radio Wired: via supplied USB cable
Data Transmission Range	Secure transmission of data between reader and specified PC, laptop or PDA up to a distance of 100 meters

Ordering Guide

Item No.	Description	Qty
3122	VisionMate Wireless 2D Barcode Reader with Bluetooth Technology (includes wireless barcode reader with CodeXML software, USB cable and Belkin Bluetooth device with associated software)	1 each

© 2008 Thermo Fisher Scientific Inc. All rights reserved.
Bluetooth is a registered trademark of Bluetooth SIG, Inc.
Microsoft and Windows are registered trademarks of Microsoft Corporation in the United States and other countries.
All other trademarks are the property of Thermo Fisher Scientific Inc. and its subsidiaries.

North America: Tel: 800.345.0206 | email: matrix.info@thermofisher.com
Europe: Tel: +44 (0) 161 486 2110 | email: matrix.eu.info@thermofisher.com
Asia: email: matrix.ap.info@thermofisher.com

web: www.thermo.com/matrix

Thermo
SCIENTIFIC